



## Anti-SRT Tag Mouse Monoclonal Antibody(3F7)

Catalog No.	EM33037-01	EM33037-02
Size	50 $\mu$ l	100 $\mu$ l

The SRT-tag is a peptide epitope comprised of 10 amino acids(TFIGAIATDT) derived from the crystalline surface layer protein of Rickettsia typhi1.

**Molecular Weight:** N/A


**Application :** WB 1:3,000-5,000 Optimal dilutions should be determined by the end user.

**Specificity:**The SRT Tag Mouse Monoclonal Antibody can detects SRT Tag recominant protein.

**Form:** Liquid ,1mg/ml

**Storage buffer:** PBS, pH 7.4, containing 0.02% sodium azide as Preservative and 50% Glycerol

**Storage:** Store at -20°C. Do not aliquot the antibody.

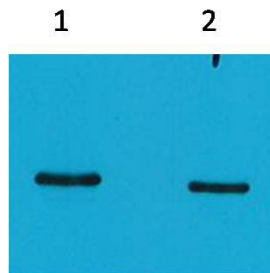

Western blot analysis of SRT-Recombinant protein with SRT mAb diluted at  
1) 1:5,000      2) 1:10,000